

**THROMBOLEX™**  
INNOVATIVE ENDOVASCULAR CATHETERS

**BASHIR™+ : A POWERFUL DEVICE TO HELP CONQUER VTE**



**BASHIR™+**  
*endovascular catheter*

Discover the family of catheters that will change your practice:  
[www.thrombolex.com](http://www.thrombolex.com)

Caution: Federal (USA) law restricts these devices to sale by or on the order of a physician (Rx Only). Prior to use, please refer to the Instructions for Use for complete product intended use and indications for use, contraindications, warnings, precautions, potential adverse events, and detailed instructions for use. Renderings for illustrative purposes. Photographs taken and on file at Thrombolex, Inc. Please contact your local Thrombolex representative for more information.

# TAP030 Rev. C

**BASHIR™**  
*endovascular catheter*

A Simplified Solution for  
Complete Thrombus Resolution



**X - HOPEFUL HEALER. ACTIVE VTE CONQUEROR.**

**THROMBOLEX™**  
INNOVATIVE ENDOVASCULAR CATHETERS

**The truth is, VTE affects more than just the patient.** Current treatments may not provide a complete solution, placing a burden on your budget and staff. Thrombolex has the solution.

The BASHIR™ endovascular family of catheters has innovative expandable infusion baskets for treatment of arterial and venous thrombotic conditions. Its hybrid design enables both mechanical and pharmacological treatment, to achieve optimal results.

Discover the family of catheters that will change your practice: [thrombolex.com](http://thrombolex.com).

## THE BASHIR™ ENDOVASCULAR CATHETER DIFFERENCE

### THE BASHIR™ ENDOVASCULAR CATHETER BASKET EXPANSION\*

#### HOW THE BASHIR™ ENDOVASCULAR CATHETER WORKS

1. Immediate reperfusion via basket expansion and the creation of cross-sectional fissures.
2. Targeted delivery of solutions, including lytics, through 6 mini-infusion catheters.
3. Blood flow restoration using endogenous and exogenous lytics to rapidly accelerate lysis.



Atraumatic nitinol infusion basket expands from 3mm up to 45mm  
7 Fr sheath recommended

\*BASHIR™ Endovascular Catheter actual size

### BASHIR™ ENDOVASCULAR CATHETER TREATMENT SCENARIOS

#### Venous thrombus

- ilio caval (IVC filter)
- iliofemoral
- femoral
- Paget-Schroetter syndrome

#### Arterial thrombus

- aortic
- iliac
- pulmonary artery (BASHIR™ N-X Endovascular Catheter)
- femoral

## BASHIR™ endovascular catheter


### CASE STUDIES




Pre-intervention venogram of occluded IVC filter





Post-intervention (BASHIR™ Endovascular Catheter+40), total 25mg r-tPA

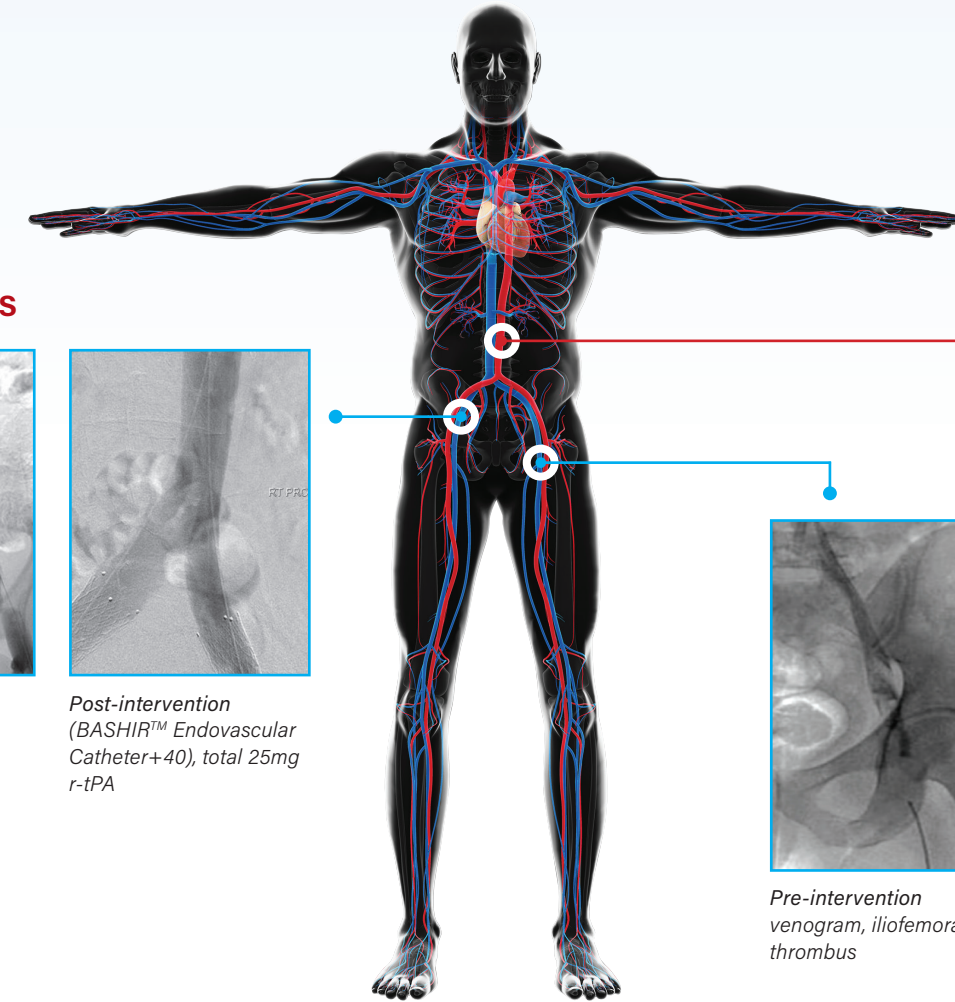


Pre-intervention aortogram

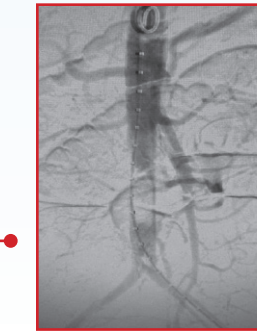

Post-intervention 30 minutes post-intervention BASHIR™ Endovascular Catheter with 2 x 5mg boluses r-tPA

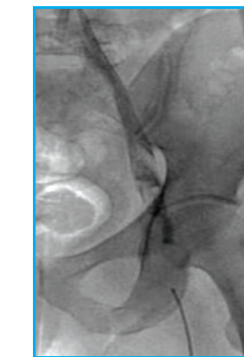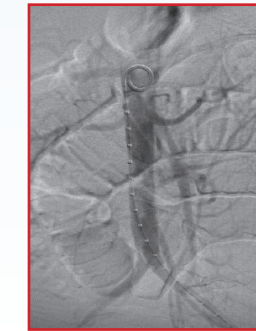

Pre-intervention venogram, iliofemoral thrombus

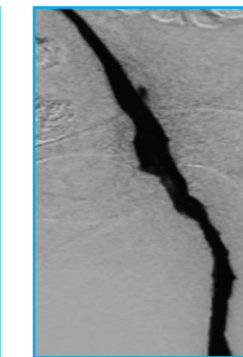

Post-intervention 16mg total r-tPA

Customer Service: 844.792.6300 Email Inquiries: [INFO@THROMBOLEX.COM](mailto:INFO@THROMBOLEX.COM) Orders: [ORDERS@THROMBOLEX.COM](mailto:ORDERS@THROMBOLEX.COM) Website: [THROMBOLEX.COM](http://THROMBOLEX.COM)